

C-390



Discs containing antibiotic C-390 for the identification of *Pseudomonas aeruginosa* by the agar diffusion method




The microorganism *Pseudomonas aeruginosa* shows a specific resistance to the antibiotic C-390 compared to the other microorganisms of *Pseudomonas* genus and the other non-fermenting bacteria

AT THE END OF THE INCUBATION PERIOD, THE INHIBITION HALOS PERMIT THE DISTINCTION:

Diffusion method on MUELLER HINTON AGAR medium

Pseudomonas aeruginosa ATCC 27853

Burkholderia cepacia ATCC 1777

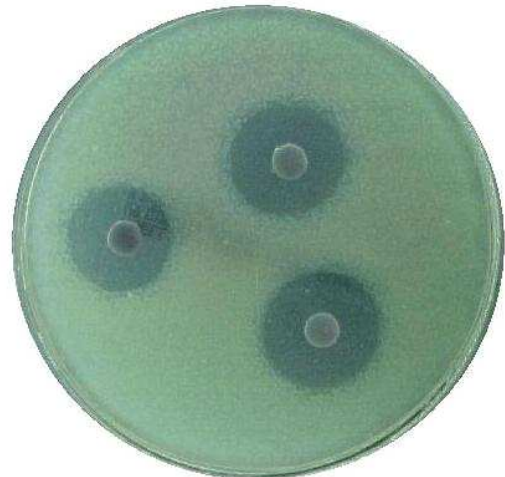

no halo:
Pseudomonas aeruginosa

halo \geq 12 mm:
Other *Pseudomonas* spp. and non fermenting



LIOFILCHEM Bacteriology Products

Via Scozia Zona Ind.le - 64026 Roseto D.A. (TE) - Italy

Phone +390858930745 Fax +390858930330

Internet: www.liofilchem.net

E-Mail: liofilchem@liofilchem.net